

MKS® Lumbago Vario

Lumbar Bandage with silicone pad stimulating blood circulation

Indications

- Severe lower back pain
- Severe lumbago
- Relapsing lumbago
- Osteochondrosis
- Spondylolysis without displacement of vertebrae
- Mild lumbar deformity (with displacement of vertebrae)
- Spondylarthrosis
- Sacroiliac joint irritation
- Lumbar muscle insufficiency


Features

- Semi-elastic bandaging material
- Lateral reinforcing stays
- Firm abdominal piece with hook and loop fastener
- Back pocket suitable for silicone pad stimulating blood circulation or MKS® Osteo back plate
- Circumference and shape can be adjusted
- Silicone pad stimulating



Mode of Operation

Due to its semi-elastic bandaging material the MKS® Lumbago Vario adjusts comfortably to the body. Through building up circular pressure the lumbar vertebral column is supported. The firm fastener element provides the necessary abdominal pressure which can be enhanced through the additional tensioning straps. The inserted silicone pad in the integrated back pocket stimulates blood circulation.

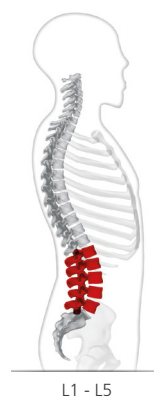

Size chart and article numbers

Size	Standard	Extra
Body circumference [cm]	70 - 95	95 - 125
Abdominal height [cm]	18	18
Back height [cm]	25	25
Art.-No.	581-069-S	581-070-S

Build-up

The MKS® Lumbago Vario consists of two side pieces (1) and (2), the middle piece (3) and the textile cover (4).

The MKS® Lumbago Vario can be adjusted for different circumferences and shapes to accommodate various physiques



Other sizes available on request

Meaning of indicated dimensions

