

CDS® Wrist Brace DUO WITH NEW HINGE

☎ 00800 03032900

Dynamic spring-loaded wrist orthosis

Patents DE 199 04 554 / EP 1 482 881 / US 7,083,583 / DE 10 2015 012 320 / EP 3 146 944

Indications

- Joint contractures due to
 - vegetative state
 - cerebral palsy
 - genetic syndroms
 - epilepsy
 - spina bifida
 - near-drowning
 - burns
 - stroke
 - spinal injury
 - paralysis
 - craniocerebral injury
 - capsular ligament injury
 - surgery
 - arthrosis and chronic polyarthritis
- to prevent recurring constracures following arthrolysis

Features

- Therapy in extension and flexion
- spring tension individually adjustable
- Stops in flexion and extension adjustable in 15° increments
- ROM range and switch point adjustable in 10° increments using 360° feature
- Locking function: Easy to put on and take off the orthosis by locking the joint in any position. Immobilisation possible
- Redression range 160°
- High wearing comfort thanks to air-permeable, lightweight aluminium shell elements and soft cover on the hand section
- High flexibility thanks to adjustable, adaptable shell elements
- Finger section can be individually adjusted to the course of treatment
- Easy handling thanks to quick-release fasteners on the forearm



Mode of operation

The CDS® Wrist Brace Duo is based on the CDS®-principle and has been designed to treat an extension and flexion deficit of the wrist. The brace applies a dynamic low load prolonged stretch to stimulate growth in the contracted tissue. The infinitely adjustable redression range prevents harmful overstretching.

Size chart and article numbers - NEW SIZES AVAILABLE

Size	L/L	L/M	L/S	M/M
Hinge	CDS® DUO Hinge			
Shells	Aluminium Shells			
Length at forearm medial [cm]*	21 - 23,5	21 - 23,5	21 - 23,5	17 - 19
Circumference at forearm [cm]	22 - 30	18 - 24	14 - 19	18 - 24
Circumference of wrist [cm]	10	8	6	8
Art.-No. left	930LL-L	930LM-L	930LS-L	930MM-L
Art.-No. right	930LL-R	930LM-R	930LS-R	930MM-R

*in combination with the respective shell set

Optional accessories

	Art.-No.
carrying strap L	580-200-S
carrying strap S (JUNIOR)	580-210-S
omosan Shoulder Brace connecting set	37-7100

Modules for the modular assembly

When combining modules of different sizes, we recommend consulting albrecht to ensure the compatibility of the modules.

Shell Sets Forearm



Size	L/L	L/M	L/S	M/M
Shells	Aluminium Shells			
Circumference at forearm [cm]	22 - 30	18 - 24	14 - 19	18 - 24
Art.-No. left	824LL-L	824LM-L	824LS-L	824MM-L
Art.-No. right	824LL-R	824LM-R	824LS-R	824MM-R

Forearm shell set XL (24-36cm): 824LXL-L / 824LXL-R
Forearm shell set XXL (28-38cm): 824LXXL-L / 824LXXL-R

Hinge CDS® Wrist Brace DUO



Hinge	Hinge Set CDS® DUO Hinge L	Hinge Set CDS® DUO Hinge M
Use for	L/L, L/M, L/S	M/M
Length at forearm medial [cm]*	21 - 23,5	17 - 19
Art.-No. left	870-215-L	870-216-L
Art.-No. right	870-215-R	870-216-R

*in combination with the respective shell set

Palm Plates



Size	L plus*	L	M	S
Width [cm]	10	10	8	6
Art.-No. left	827L-PL	827L-L	827M-L	827S-L
Art.-No. right	827L-PR	827L-R	827M-R	827S-R

* Like palm plate L, only with wide finger plating.

Maßblatt obere Extremität

Bei Bestellungen bitte kopieren und ausgefüllt senden.
E-Mail: info@albrechtgmbh.com
Fax: 08051-961 29 30

Bestellung Kostenvoranschlag Datum: _____

Sanitätshaus: _____

Kundennummer: _____

Orthopädietechniker: _____

Patienten ID: _____

Seite

Links Rechts

Produktname

- Elbo-track
- CDS® Ellbogenschiene
- CDS® Ellbogenschiene DUO
- CDS® Sup-Prosan
- CDS® Handgelenkschiene
- CDS® Handgelenkschiene DUO
- CDS® Manusan
- Manusan
- CDS® Kombinationsorthese Ellbogen-Hand
- sonstige: _____

Wirkrichtung (bei CDS®)


Ellbogen Extension Flexion

Hand Extension Flexion

Für Kinder

Farbe

- Grün 
- Orange 
- Blau 
- Pink 

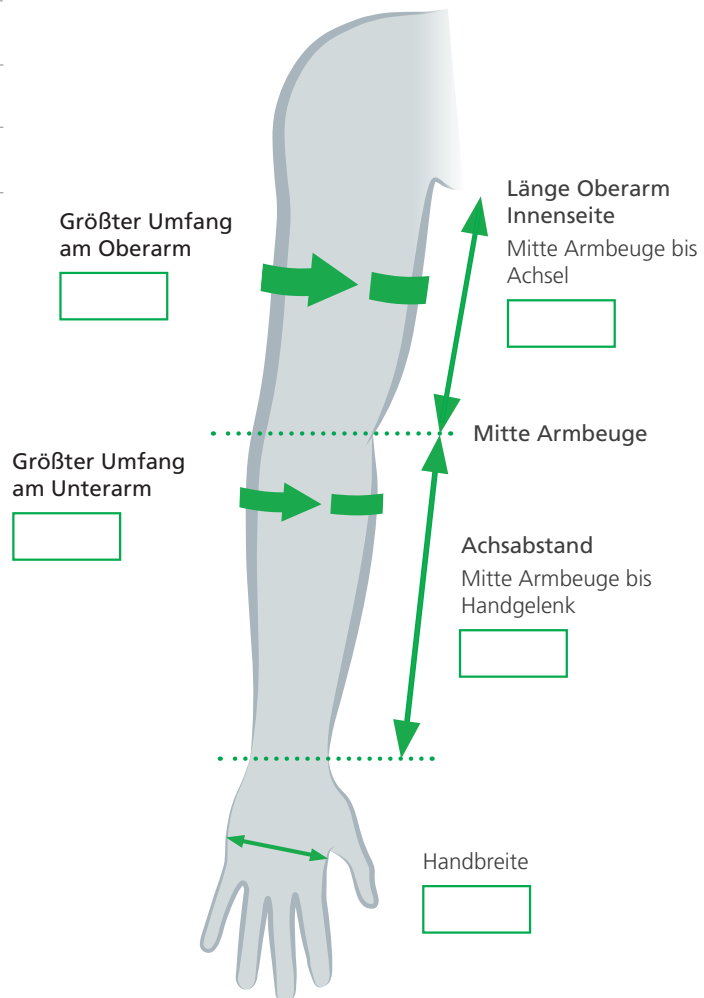

Technische Umsetzung (von albrecht ausgefüllt)

- Größenmodifikation
- modulare Änderung

Art.-Nr.: _____

Freigabe: _____