

MKS® Lumbago plus

☎ 00800 03032900

Stabilising orthosis with tensioning elements and attachment and removal options and pad

Indications


- Severe lower back pain
- Severe lumbago
- Relapsing lumbago
- Osteochondrosis
- Mild lumbar sciatica
- Mild radicular syndrome
- Spondylolysis without displacement of vertebrae
- Mild lumbar deformity (with displacement of vertebrae)
- Spondylarthrosis
- Sacroiliac joint irritation
- Lumbar muscle insufficiency

Features

- Additional tensioning straps
- Silicone Pad
- Semi-elastic bandaging material
- Lateral reinforcing stays
- Firm abdominal piece with micro hook and loop fastener
- Back pocket suitable for silicone pad stimulating blood circulation or MKS® Osteo back plate



Silicone Pad

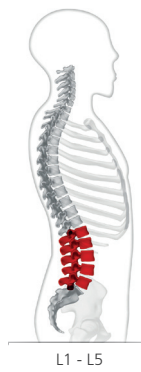

L1 - L5



Mode of operation

Due to its semi-elastic bandaging material the MKS® Lumbago plus adjusts comfortably to the body. Through building up circular pressure the lumbar vertebral column is supported. The firm fastener element provides the necessary abdominal pressure which can be enhanced through the additional tensioning straps. The inserted silicone pad in the integrated back pocket stimulates blood circulation.

Size chart and article numbers

| Size | 1 | 2 | 3 | 4 | 5 | 6 |
|-------------------------|-----------|-----------|-----------|-----------|-----------|-----------|
| Body circumference [cm] | 60 - 70 | 70 - 80 | 80 - 90 | 90 - 105 | 105 - 120 | 120 - 135 |
| Abdominal height [cm] | 18 | 18 | 18 | 18 | 18 | 18 |
| Back height [cm] | 22 | 23 | 23 | 24 | 24 | 24 |
| Art.-No. | 582-501-S | 582-502-S | 582-503-S | 582-504-S | 582-505-S | 582-506-S |

Dimension sheet torso

For orders please copy and send filled in form to:
 E-Mail: info@albrechtgmbh.com
 Fax: 08051-961 29 30

Order Quotation Date: _____

Medical supply store: _____

Customer ID: _____

Orthopaedic technician: _____

Patient ID: _____

Product name

- MKS® Body belt
- MKS® Body belt plus
- MKS® Vario
- MKS® Lumbago
- MKS® Lumbago plus
- MKS® Lumbago Vario
- MKS® Pontsana®
- MKS® Pontsana® plus
- MKS® Pontsana® mobil
- MKS® Pontsana® stabil
- MKS® Osteo
- MKS® Osteo plus
- MKS® Osteo flex
- MKS® Osteo flex plus
- MKS® Lumbo
- MKS® Lumbo forte
- MKS® Lumbo plus
- MKS® Lumbo forte plus
- MKS® Lumbo mobil
- MKS® Lumbo fenestra
- MKS® Extensa
- MKS® Extensa plus
- MKS® Thorecta
- MKS® Thorecta plus
- MKS® Core Compression Bandage
- other: _____

Abdominal height
 Maximum front height of the textile body belt

Front height
 Measured from the upper edge of the pubic bone to the sternum (3-4 below the collarbone)

Body circumference
 Measured at the middle of the torso (just above the navel)

Back height reclination

Measured from the sacrum (beginning of gluteal fold) to circa 2-3 cm below cervical vertebra C7

Back height

Measured from the sacrum to circa 2-3 cm below the tips of the shoulder blades



Technical implementation (to be filled in by albrecht)

- Size modification
- Modular adaptation

Art.-No.: _____

Approval: _____