

Orthosis for active relief and correction of the lumbar/thoracic spine in the sagittal plane with up and down adjustment option

Sizes


 This product is also available in pediatric sizes.

Indications

- Stable, osteoporotic and traumatic vertebral fractures
- Pain caused by poor posture, e.g. rounded shoulders
- Spondylitis
- Tumor
- Condition after intervertebral disc surgery
- Muscle insufficiency
- Spondylodiscitis
- Abdominal muscle insufficiency
- Sintered vertebral fractures in osteoporosis
- Osteoporosis
- Poor posture
- Hollow round back
- Chronic lumbago caused by hyperlordosis
- Facet syndrome
- Degenerative changes in the lumbar spine

Features

- Semi-elastic bandaging material
- Lateral reinforcing staysäbe
- Firm abdominal piece with micro hook and loop fastener
- Adjustable back plate with incorporated aluminium profile
- Individual strap course
- Reclination straps

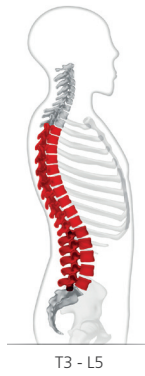

Mode of operation

The MKS® Osteo functions according to the reclination principle.

The two shoulder straps attached to the back plate straighten the vertebral column to relieve the thoracic spine. The wide lumbar bandage straightens and hence relieves the lumbar vertebral column. The reclination straps are attached to the lumbar bandage through additional lateral straps. Hereby rotation is limited and stabilisation is ensured.

Size chart and article numbers

| Size | | 1 | 2 | 3 | 4 | 5 | 6 |
|-------------------------------------|------|---------|---------|---------|----------|-----------|-----------|
| Body circumference [cm] | | 60 - 70 | 70 - 80 | 80 - 90 | 90 - 105 | 105 - 120 | 120 - 135 |
| Abdominal height [cm] | 18 | | | | | | |
| Back height reclination [cm] | S 46 | 782C-1 | 782C-2 | 782C-3 | 782C-4 | 782C-5 | 782C-6 |
| | M 50 | 782A-1 | 782A-2 | 782A-3 | 782A-4 | 782A-5 | 782A-6 |
| | L 56 | 782B-1 | 782B-2 | 782B-3 | 782B-4 | 782B-5 | 782B-6 |

Reduced height

| | | | | | | | |
|-------------------------------------|-------|---------|---------|---------|---------|---------|---------|
| Abdominal height [cm] | 14 | | | | | | |
| Back height reclination [cm] | XS 43 | 782D-1N | 782D-2N | 782D-3N | 782D-4N | 782D-5N | 782D-6N |
| | S 46 | 782C-1N | 782C-2N | 782C-3N | 782C-4N | 782C-5N | 782C-6N |
| | M 50 | 782A-1N | 782A-2N | 782A-3N | 782A-4N | 782A-5N | 782A-6N |
| | L 56 | 782B-1N | 782B-2N | 782B-3N | 782B-4N | 782B-5N | 782B-6N |

Modules for the modular assembly

| Name | Art.-No. | Included as standard in the products Art.-No. 782... | | | | | | | | | | | | | | | | | |
|---|-----------|--|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|
| | | C1 | A1 | B1 | C2 | A2 | B2 | C3 | A3 | B3 | C4 | A4 | B4 | C5 | A5 | B5 | C6 | A6 | B6 |
| Dorsal plate with reclination straps | | | | | | | | | | | | | | | | | | | |
| S 46 cm | 761-066-S | • | | | • | | | • | | | • | | | • | | | • | | |
| M 50 cm | 761-067-S | | • | | | • | | | • | | | • | | | • | | | • | |
| L 56 cm | 761-068-S | | | • | | | • | | | • | | | • | | | • | | | • |
| Body belt Gr. 1 60 - 70 cm | 582-001-S | • | • | • | | | | | | | | | | | | | | | |
| Body belt Gr. 2 70 - 80 cm | 582-002-S | | | | • | • | • | | | | | | | | | | | | |
| Body belt Gr. 3 80 - 90 cm | 582-003-S | | | | | | | • | • | • | | | | | | | | | |
| Body belt Gr. 4 90 - 105 cm | 582-004-S | | | | | | | | | | • | • | • | | | | | | |
| Body belt Gr. 5 105 - 120 cm | 582-005-S | | | | | | | | | | | | | • | • | • | | | |
| Body belt Gr. 6 120 - 135 cm | 582-006-S | | | | | | | | | | | | | | | | • | • | • |



Reduced height

| Name | Art.-No. | Included as standard in the products Art.-No. 782... | | | | | | | | | | | | | | | | | | | | | | | | |
|---|------------|--|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|----|---|
| | | D1 | C1 | A1 | B2 | D2 | C2 | A2 | B2 | D3 | C3 | A3 | B3 | D4 | C4 | A4 | B4 | D5 | C5 | A5 | B5 | D6 | C6 | A6 | B6 | |
| Dorsal plate with reclination straps | | | | | | | | | | | | | | | | | | | | | | | | | | |
| XS 43 cm | 761-065-S | • | | | | • | | | | • | | | | • | | | | • | | | • | | | | • | |
| S 46 cm | 761-066-S | | • | | | | • | | | | • | | | | • | | | | • | | | • | | | • | |
| M 50 cm | 761-067-S | | | • | | | | • | | | | • | | | | • | | | | • | | | | • | | |
| L 56 cm | 761-068-S | | | | • | | | | • | | | | • | | | | • | | | | • | | | | • | |
| Body belt Gr. 1 60 - 70 cm | 582-001N-S | • | • | • | • | | | | | | | | | | | | | | | | | | | | | |
| Body belt Gr. 2 70 - 80 cm | 582-002N-S | | | | | • | • | • | • | | | | | | | | | | | | | | | | | |
| Body belt Gr. 3 80 - 90 cm | 582-003N-S | | | | | | | | | • | • | • | • | | | | | | | | | | | | | |
| Body belt Gr. 4 90 - 105 cm | 582-004N-S | | | | | | | | | | | | | • | • | • | • | | | | | | | | | |
| Body belt Gr. 5 105 - 120 cm | 582-005N-S | | | | | | | | | | | | | | | | | | • | • | • | • | | | | |
| Body belt Gr. 6 120 - 135 cm | 582-006N-S | | | | | | | | | | | | | | | | | | | | | | • | • | • | • |



Dimension sheet torso

For orders please copy and send filled in form to:
 E-Mail: info@albrechtgmbh.com
 Fax: 08051-961 29 30

Order Quotation Date: _____

Medical supply store: _____

Customer ID: _____

Orthopaedic technician: _____

Patient ID: _____

Product name

- MKS® Body belt
- MKS® Body belt plus
- MKS® Vario
- MKS® Lumbago
- MKS® Lumbago plus
- MKS® Lumbago Vario
- MKS® Pontsana®
- MKS® Pontsana® plus
- MKS® Pontsana® mobil
- MKS® Pontsana® stabil
- MKS® Osteo
- MKS® Osteo plus
- MKS® Osteo flex
- MKS® Osteo flex plus
- MKS® Lumbo
- MKS® Lumbo forte
- MKS® Lumbo plus
- MKS® Lumbo forte plus
- MKS® Lumbo mobil
- MKS® Lumbo fenestra
- MKS® Extensa
- MKS® Extensa plus
- MKS® Thorecta
- MKS® Thorecta plus
- MKS® Core Compression Bandage
- other: _____

Abdominal height
 Maximum front height of the textile body belt

Front height
 Measured from the upper edge of the pubic bone to the sternum (3-4 below the collarbone)

Body circumference
 Measured at the middle of the torso (just above the navel)



Back height reclination

Measured from the sacrum (beginning of gluteal fold) to circa 2-3 cm below cervical vertebra C7

Back height

Measured from the sacrum to circa 2-3 cm below the tips of the shoulder blades

Technical implementation (to be filled in by albrecht)

- Size modification
- Modular adaptation

Art.-No.: _____

Approval: _____