

omosan Shoulder Brace

Shoulder joint orthosis for physiological guidance and stabilisation

Indications

- Subluxation of the shoulder in hemiplegia
- Paralysis of the muscles surrounding the shoulder
- Lesion of the brachial plexus
- Peripheral nerve damage in the cervical spine and shoulder region
- Mild to moderate instability of the shoulder joint
- Conservative treatment for habitual shoulder dislocation
- Inflammation
- Degeneration

Features

- Partially elastic bandage material
- Shoulder-encompassing element
- Forearm cuff with adjustable straps
- One size
- Can be used on the right or left side
- Can be combined with CDS® braces for the upper extremities using the **omosan Shoulder Brace** connecting set.



Mode of operation

The **omosan Shoulder Brace** adapts well to the body thanks to its partially elastic bandage material. The material surrounding the shoulder, which is fastened on the other side via the armpit with a strap, ensures that the humeral head is returned to its original position and prevents subluxation. The forearm cuff with its adjustable straps ensures that the entire arm is pulled towards the head.

Size chart and article numbers

Name	omosan Shoulder Brace	omosan Shoulder Brace connecting set
Size	One size fits all	One size fits all
Art.-No.:	37-7000	37-7100

The respective drilling templates can be found here



omosan Shoulder Brace connecting set

omosan Shoulder Brace offers new (connection-) possibilities

With the omosan shoulder brace in combination with the CDS® orthoses for the upper extremities, patients can be treated individually from the fingers to the shoulder joint.

CDS® Wrist Brace
Extension / Flexion



CDS® Elbow Brace
Extension / Flexion



CDS® Wrist Brace DUO
Extension / Flexion



CDS® Elbow Brace DUO
Extension / Flexion



CDS® Manusan



CDS® Combination Brace
Elbow-wrist
Extension / Flexion



Dimension sheet upper extremity

For orders please copy and send filled in form to:
 E-Mail: info@albrechtgmbh.com
 Fax: 08051-961 29 30

Order Quotation Date: _____

Medical supply store: _____

Customer ID: _____

Orthopaedic technician: _____

Patient ID: _____

Side

Left Right

Product name

- Elbo-track
- CDS® Elbow Brace
- CDS® Elbow Brace DUO
- CDS® Sup-Prosan
- CDS® Wrist Brace
- CDS® Wrist Brace DUO
- CDS® Manusan
- Manusan
- CDS® Combination Brace Elbow-Wrist
- other: _____

Direction of spring force (for CDS®)


Elbow Extension Flexion

Hand Extension Flexion

For children

Colour

- Green 
- Orange 
- Blue 
- Pink 

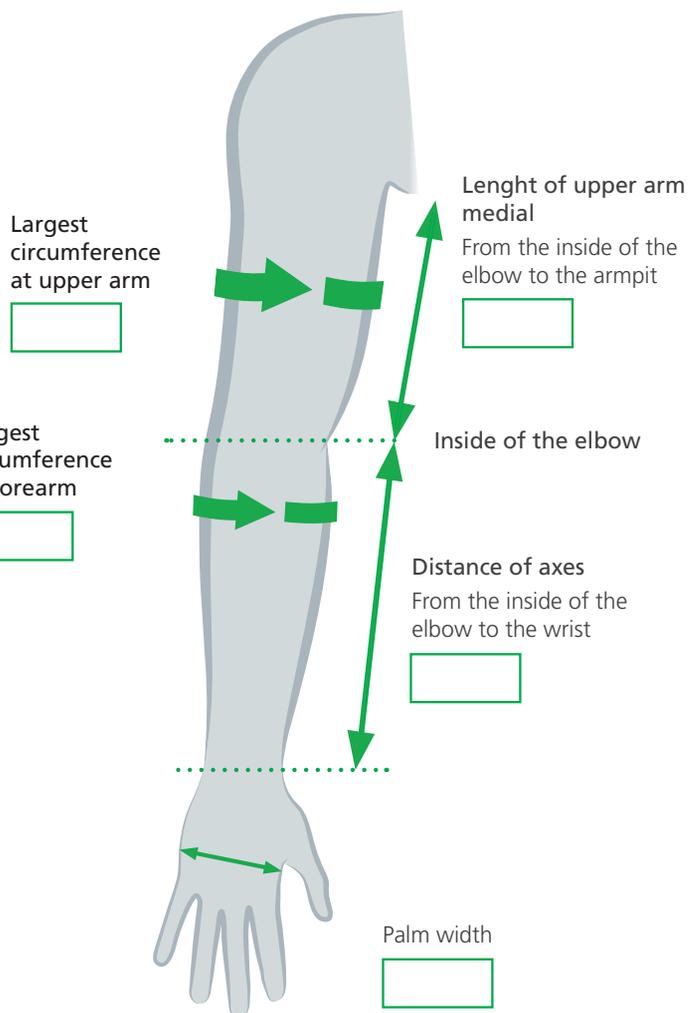

Technical implementation (to be filled in by albrecht)

- Size modification
- Modular adaptation

Art.-No.: _____

Approval: _____